

GRADIA CHARACTERIZATION

SHOULDER DENTIN						
	SD2	SD3	SD4	SD5	SD7	SD8
OPAQUS DENTIN INTENSIVE						
	ODI1	ODI2	ODI3	ODI4	ODI5	ODI6
ENAMEL INTENSIVE						
	EI0	EI1	EI3	EI5		
HALO ENAMEL						
	HE1					
PEARL ENAMEL						
	PE1	PE3				
TRANSLUCENT						
	T0	T1	T2	T4	T5	
CERVICAL TRANSLUCENT						
	CT2	CT3	CT4			
MAMELON STAIN						
	MS2	MS3	MS5			
INTENSIVE COLOR						
	IC0	IC1	IC2	IC3	IC4	IC5
FLOWABLE COMPOSITE						
	F-WT	F-WE	F-T0			

IRRADIATION TIME CHART


	LABOLIGHT LV-III, II		STEPLIGHT SL-I
	Pre-cure	Final Cure	Pre-cure
FOUNDATION OPAQUE OPaque MARGIN OPAQUE GUM OPAQUE GUM OPAQUE MODIFIER	1 min.	—	—
DENTIN OPAQUS DENTIN OPAQUS DENTIN INTENSIVE SHOULDER DENTIN ENAMEL ENAMEL INTENSIVE HALO ENAMEL PEARL ENAMEL TRANSLUCENT CERVICAL TRANSLUCENT GUM GUM MODIFIER GUM TRANSLUCENT	30 sec.	3 min.	10 sec.*
INTENSIVE COLOR MAMELON STAIN FLOWABLE COMPOSITE	30 sec.	—	10 sec.*

*For one surface of a single crown


GRADIA GUM SHADE COMBINATION CHART

OPAQUE	GUM OPAQUE	 GO11	 GO12	 GO13					GUM FIBER	 GF71
	GUM OPAQUE MODIFIER	 GOM51								
BODY (Paste)	GUM	 G20	 G21	 G22	 G23	 G24				
MODIFIER (Gel)	GUM MODIFIER	 GM30	 GM31	 GM32	 GM33	 GM34	 GM35	 GM36		
TRANSLUCENT (Gel)	GUM TRANSLUCENT	 GT41								

Example of Build-up

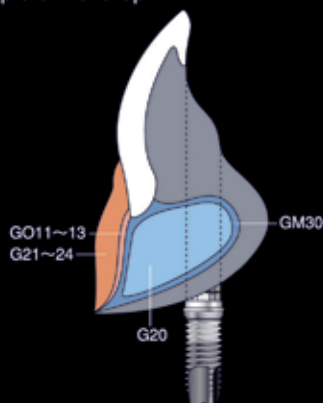

Example of Multiple-Layering Technique

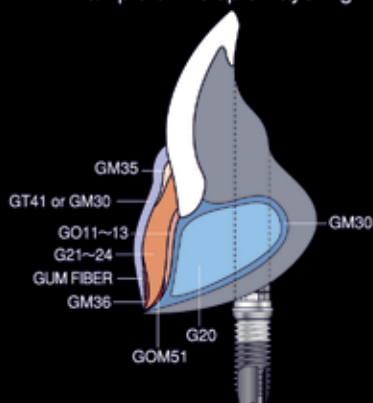

GC EUROPE N.V.
Head Office
Researchpark Haasrode-Leuven 1240
Interleuvenlaan 33
B - 3001 Leuven
Tel. +32.16.74.10.00
Fax. +32.16.40.48.32
info@gceurope.com
<http://www.gceurope.com>

GC

GRADIA

GRADIA

LIGHT-CURED COMPOSITE FOR CROWN & BRIDGE,
INLAYS AND VENEERS

GRADIA

GUM SHADES LIGHT CURED GUM SHADE COMPOSITE SYSTEM

SHADE COMBINATION CHART

		A1	A2	A3	A3.5	A4	B1	B2	B3	B4	C1	C2	C3	C4	D2	D3	D4	A0	B0	
OPAQUE	FOUNDATION OPAQUE	 FO																 FO		
	MARGIN OPAQUE	 MO																—		
	OPAQUE ORIGINAL SHADE	 OA1	 OA2	 OA3	 OA3.5	 OA4													 OW	
	OPAQUE CLASSICAL SHADE	 OA1-C	 OA2-C	 OA3-C	 OA3.5-C	 OA4-C		 OB1	 OB2	 OB3	 OB4	 OC1	 OC2	 OC3	 OC4	 OD2	 OD3	 OD4		
CERVICAL	OPAQUE DENTIN, OPAQUE DENTIN INTENSIVE*1 ORIGINAL SHADE	 ODA2	 ODA3	 ODA3.5	 ODA4	 ODI3													 ODA1	 ODB1
	OPAQUE DENTIN, OPAQUE DENTIN INTENSIVE*1 CLASSICAL SHADE	 ODA2-C	 ODA3-C	 ODA3.5-C	 ODA4-C		 ODB2	 ODB3	 ODB4	 ODI2	 ODC2	 ODC3	 ODC4	 ODI5	 ODD3	 ODD4	 ODI6			
DENTIN	OPAQUE DENTIN ORIGINAL SHADE	 ODA1	 ODA2	 ODA3	 ODA3.5	 ODA4													 ODA0	 ODB0
	OPAQUE DENTIN CLASSICAL SHADE	 ODA1-C	 ODA2-C	 ODA3-C	 ODA3.5-C	 ODA4-C	 ODB1	 ODB2	 ODB3	 ODB4	 ODC1	 ODC2	 ODC3	 ODC4	 ODD2	 ODD3	 ODD4			
	DENTIN ORIGINAL SHADE	 DA1	 DA2	 DA3	 DA3.5	 DA4													 DA0	 DB0
	DENTIN CLASSICAL SHADE	 DA1-C	 DA2-C	 DA3-C	 DA3.5-C	 DA4-C	 DB1	 DB2	 DB3	 DB4	 DC1	 DC2	 DC3	 DC4	 DD2	 DD3	 DD4			
ENAMEL	ENAMEL	 E2		 E3		 E4	 E1	 E2		 E3	 E2	 E3	 E4	 E2	 E3			 E0		
TRANSLUCENT	ENAMEL INTENSIVE	 EI1																 EI0		
	CERVICAL TRANSLUCENT	*2	 CT2	 CT4	 CT3		*2	 CT2	 CT4		*2	 CT2	 CT3	 CT2	 CT4			—		

*1: OPAQUE DENTIN and OPAQUE DENTIN INTENSIVE are used as cervical colors.

*2: The use of ENAMEL INTENSIVE (EI1) will give almost the same effects as CERVICAL TRANSLUCENT.

